

FIXTER

STERNAL CLOSURE PLATE

Barcode NO	Product CODE	SIZE	Label NAME	RANGE
8682091000408	SCS 2012	20x12mm	20x12mm Sternum Closure System	20-25mm
8682091000415	SCS 2512	25x12mm	25x12mm Sternum Closure System	25-30mm
8682091000422	SCS 3012	30x12mm	30x12mm Sternum Closure System	30-35mm
8682091000439	SCS 3512	35x12mm	35x12mm Sternum Closure System	35-40mm
8682091000446	SCS 4012	40x12mm	40x12mm Sternum Closure System	40-45mm
8682091000453	SCS 4512	45x12mm	45x12mm Sternum Closure System	45-50mm
8682091000460	SCS 5012	50x12mm	50x12mm Sternum Closure System	50-55mm



AN INNOVATION
FOR **CARDIO THORACIC**
AND **VASCULAR SURGERY**





Yayla-Med Tibbi Cihazlar Medikal Urunler Tic. Ltd. Sti.
Mansuroglu Mah. Ankara Cad. No:47 K:4 D:410 Bayrakli / Izmir / Turkey Zip Code: 35535

Tel: +90 232 339 88 98
www.yaylamed.com.tr

Fax: +90 232 339 88 97
info@yaylamed.com.tr

FIXTER

STERNAL CLOSURE PLATE

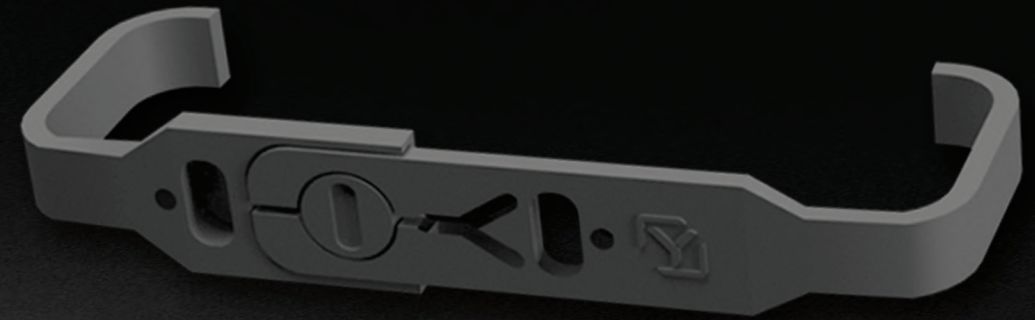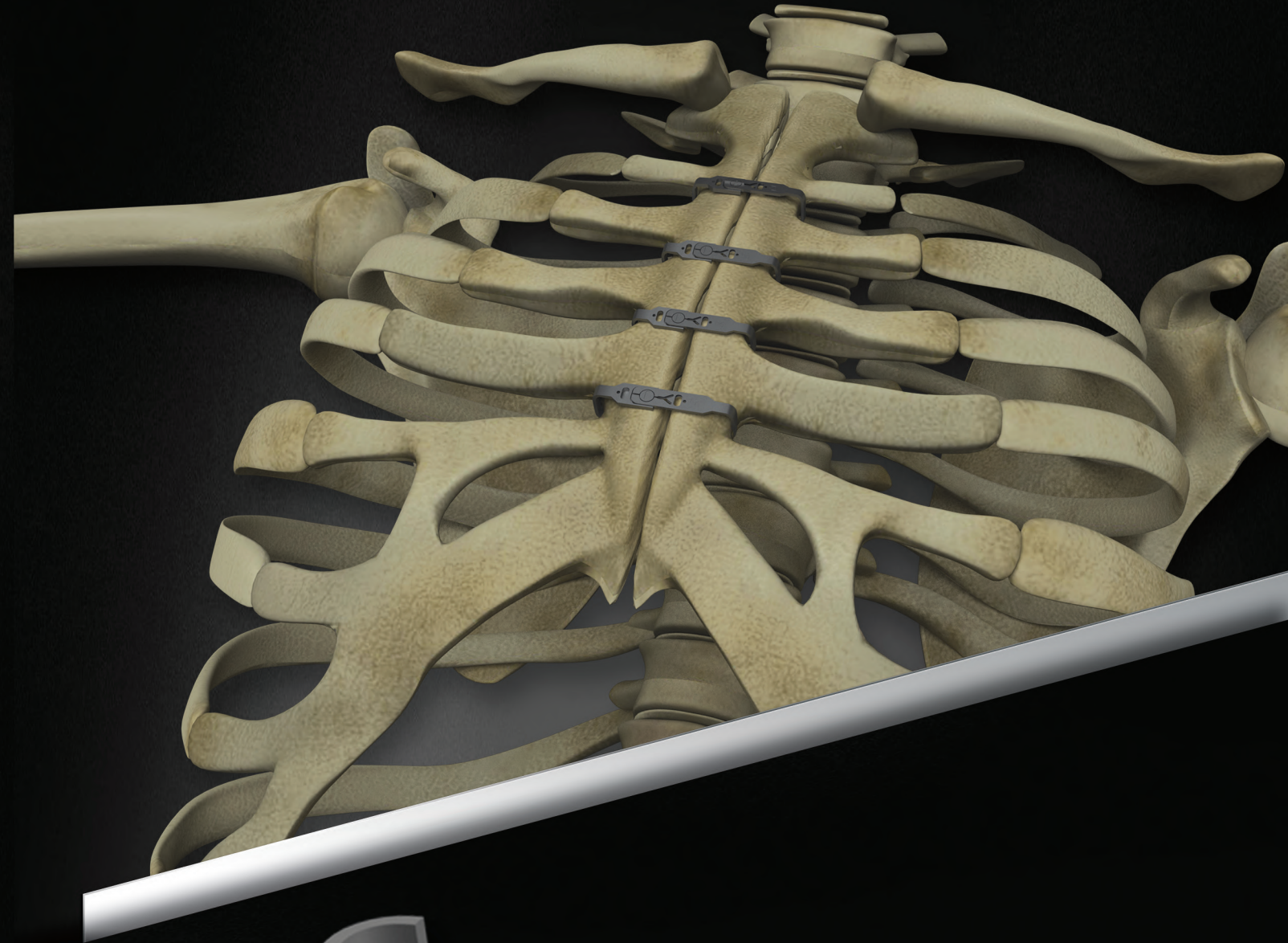

AN INNOVATION
FOR **CARDIO THORACIC**
AND **VASCULAR SURGERY**



AN INNOVATION FOR CARDIO THORACIC AND VASCULAR SURGERY

FIXTER allows you to close the sternum safely and securely. With this product, two halves of the split sternum are pulled together for rigid fixation while bone penetration is prevented at the same time.

The special design of the system ensures that the forces involved are distributed across a large bone surface – a significant advantage, especially with regard to high-risk patients. The revolutionary **FIXTER** is made of biocompatible titanium. It is easy to implant and can be opened rapidly and easily in a case of emergency immediately after the operation. What facilitates the surgeon's work a great deal on the one hand, means higher comfort for the patients on the other hand.

A firm sternal closure, stable under load immediately after the intervention, less complications and usually a faster return to normal life – these benefits translate into a real plus in the quality of life for any patient.

EASY APPLICATION AND PATIENT-FRIENDLY





The **FIXTER** is available in one basic version. **FIXTER** consist of two parts that are simply pushed together to close the sternum. Intraoperative placement involves just three simple steps. Thus, the sternum can be securely closed in a fast and uncomplicated process.



WITH JUST A TURN OPENING AND CLOSING OF THE IMPLANTS


A special ratchet mechanism ensures that the **FIXTER** implants can be opened and closed with ease. Depending on the position of the two-position screw on the back of the implant, the mechanism engages or unlocks without extra control.

HIGH SAFETY FOR HIGH RISK PATIENTS

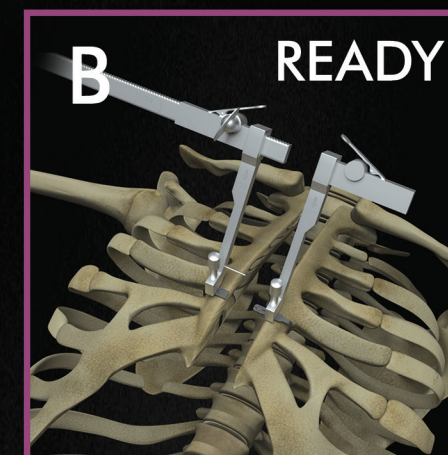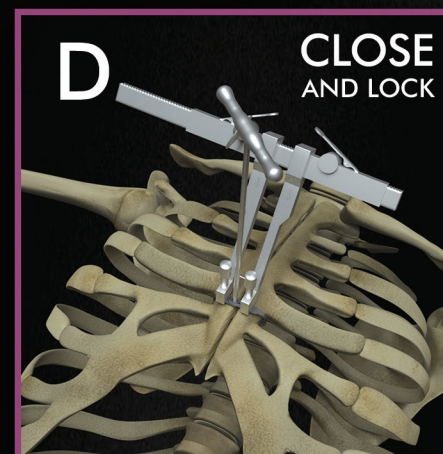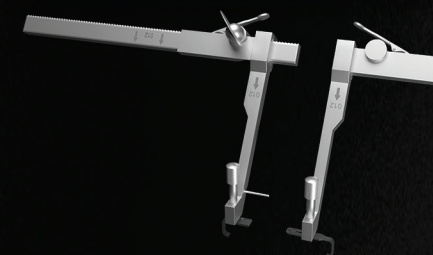
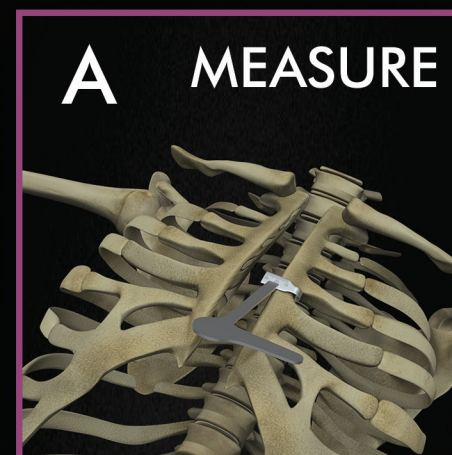


Width of the sternum measured exactly to determine the correct size of the implants to be selected.

The two-position screw of the **FIXTER** is rotated to the Closing position. This engages the ratchet mechanism. A reduction clamp is then used to press the two halves together. Once the implant has been firmly closed, the two-position screw is turned to the Locked position. This locks the implant securely.

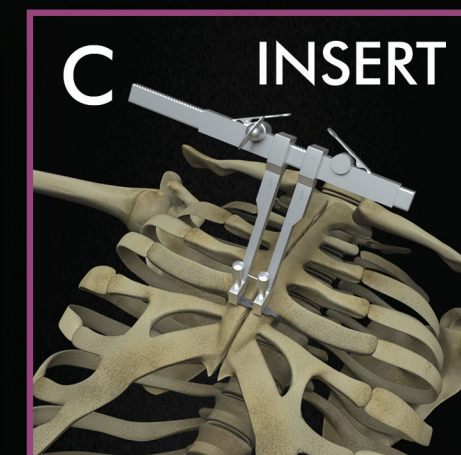

The opened **FIXTER** is inserted and then aligned parallel to the sternum.

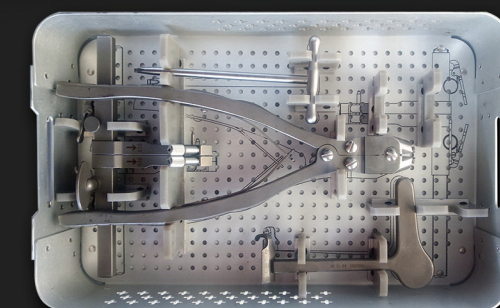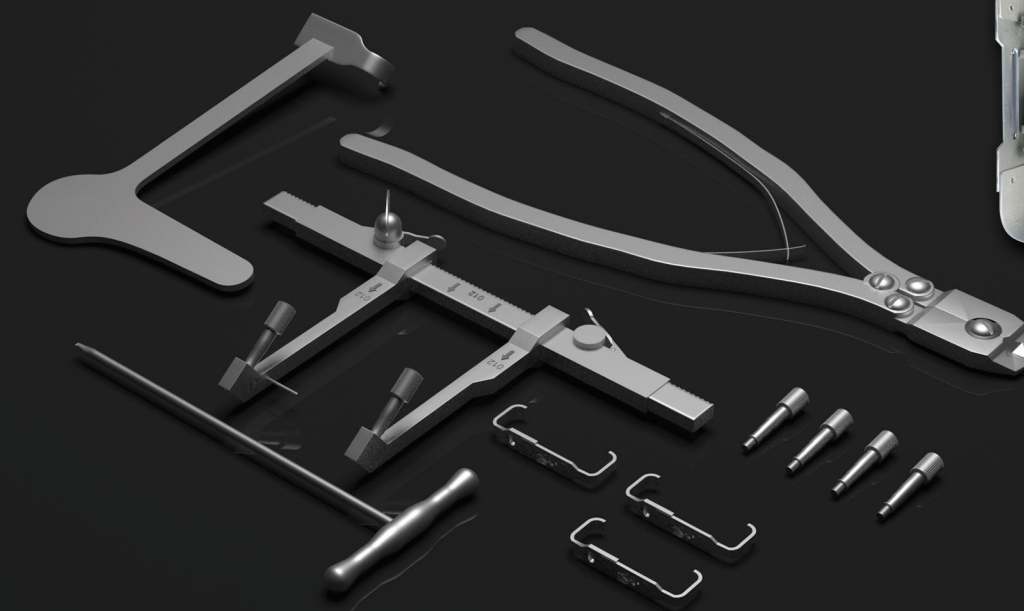

Caliper	2 Pcs
Screwdriver	1 Pcs
Cutter	1 Pcs
Tool Box	1 Pcs

FIXTER
STERNAL CLOSURE PLATE